

A simple, robust, and low-cost method to produce the PURE cell - free system.

Barbora Lavickova and Sebastian J. Maerkl*

Institute of Bioengineering,

School of Engineering

École Polytechnique Fédérale de Lausanne

Lausanne, Switzerland

Abstract

We demonstrate a simple, robust, and low-cost method for producing the PURE cell-free transcription-translation system. Our OnePot PURE system achieved a protein synthesis yield of $156 \mu\text{g}/\text{mL}$ at a cost of $0.09 \text{ USD}/\mu\text{L}$, leading to a 14-fold improvement in cost normalized protein synthesis yield over existing PURE systems. The OnePot method makes the PURE system easy to generate and allows it to be readily optimized and modified.

* Correspondence: sebastian.maerkl@epfl.ch

INTRODUCTION

Cell-free transcription-translation systems have become popular for molecular engineering [1–6]. Cell-free systems can be categorized into two main classes: cell extract and recombinant systems. Cell extracts are highly functional but complex and undefined cell-free systems. In 2001, Shimizu *et al.* demonstrated that a defined cell-free system called the “PURE” system (protein synthesis using recombinant elements) could be reconstituted from purified recombinant components [7]. Because of its defined and minimal nature, PURE is an appealing choice for biological systems engineering. The PURE system has been used for genetic network engineering [2], recombinant DNA replication [8], molecular diagnostics [9], therapeutics [10], and educational kits [11]. The PURE system also represents a viable starting point for generation of an artificial cell [12, 13] and its composition has been optimized [14, 15] and extended [16] to achieve higher functionality.

Unfortunately, producing PURE is an arduous and costly process, requiring 36 individual medium to large-scale protein purifications. PURE is now commercially available (PURExpress, New England Biolabs (NEB)), but the high-cost of the commercial system at 1.36 USD/ μ L still limits its use. Although NEB provides a few different formulations of the PURE system, the commercial system can't be customized or optimized by the user, and the precise formulation of the commercial PURE system is not publicly available. It was recently demonstrated that the PURE system could be produced using synthetic microbial “consortia” (TraMOS PURE) [17], which simplified the process of making PURE by co-expressing multiple protein components in a single *E. coli* clone combined with co-culturing of multiple strains. TraMOS PURE achieved only a ~20% protein yield compared to the commercial PURExpress and production cost was only marginally reduced from 1.36 USD/ μ L to 0.96 USD/ μ L. An earlier approach used MAGE to His-tag most PURE protein components in their endogenous locus and co-purified them from 6 strains to generate an ensemble PURE system (ePURE) [18]. The approach led initially to only minimal protein synthesis activity, and an optimized ePURE system ultimately reached a 11% protein yield compared to the original PURE system [7].

Here we present a simple, robust, and low-cost method for producing the PURE system. Our method co-cultures and induces all 36 protein producing *E. coli* clones in a single flask followed by a single Ni-NTA purification. Our “OnePot” method produces PURE at a cost

of 0.09 USD/ μ L and a protein synthesis capacity of 156 μ g/mL, which is as high as the commercial PURE system. A single batch prepares enough proteins for a total of 15 mL of PURE which is sufficient material for \sim 1,500 10 μ L reactions and can be generated in 3 days. The method produces consistent PURE across different batches and allows the optimization of individual protein components.

RESULTS AND DISCUSSION

The PURE system consists of several different components [7], that can be separated into three main categories: proteins (transcription, translation, and energy regeneration), ribosomes, and small molecule components (salts, buffers, NTPs, creatine phosphate, and folic acid). In this work, we developed a “OnePot” method for the preparation of all 36 protein components using a single mixed co-culture and Ni-NTA affinity purification step to simplify the process and decrease the cost of the PURE system. All 36 *E. coli* expression clones are cultured individually in small volumes overnight, which are then combined to inoculate a single 500 mL culture. The mixed culture is allowed to outgrow and is induced, followed by pelleting, lysis and loading of the lysate onto a Ni-NTA column for protein purification. To keep the final cost of the PURE system as low as possible, we also prepared ribosome and energy solutions (Fig. 1a, Supplementary Fig. S1). The entire process of OnePot PURE system preparation, including protein and ribosome purification and energy solution preparation, requires 4 days with 20 hours of hands-on time (Supplementary Table S1, S2, S3). To date no method has been presented in which all PURE proteins were prepared using a single co-culture and purification step [17, 18]. Moreover, other simplified protocols resulted in low protein synthesis activity as compared to the original PURE system [17].

We explored whether it is possible to adjust the protein component ratios in the OnePot PURE system simply by varying the ratios of the inoculation culture volumes added to the mixed co-culture (Supplementary Table S4). Besides ribosomal proteins, elongation factor thermo unstable (EF-Tu) is one of the most abundant proteins in rapidly growing *E. coli* [19] and it was shown to be one of the key factors for *in vitro* protein synthesis [14]. When preparing the PURE system based on the Shimizu protocol, which we refer to as “HomeMade” PURE, we observed that small concentration differences in many of

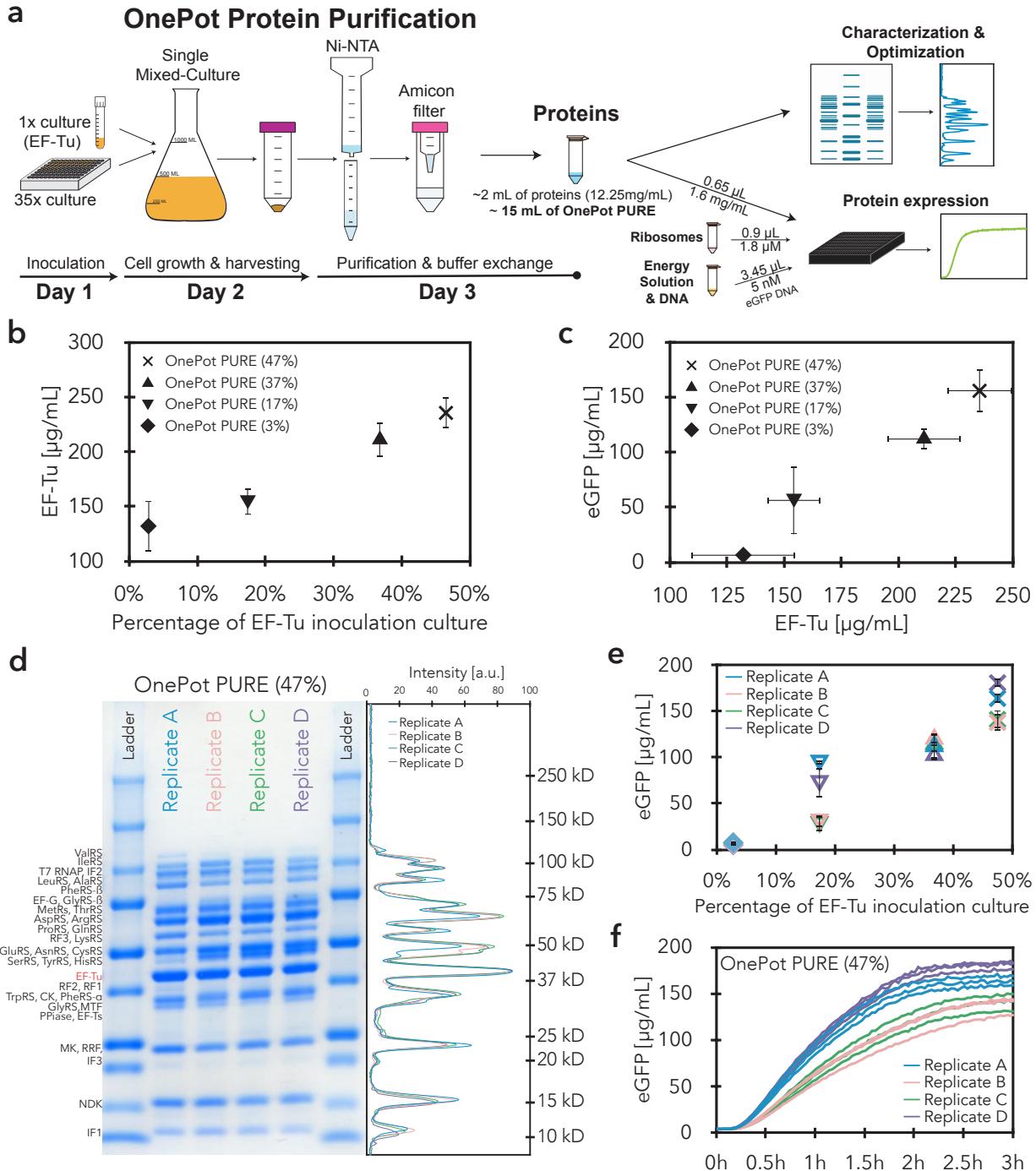


Figure 1: Caption on next page

Figure 1: OnePot PURE preparation and optimization. **(a)** All 36 PURE protein components were produced using the OnePot method, which consists of a single co-culture and a single Ni-NTA affinity purification. Different OnePot systems were produced by varying the ratio of inoculation culture EF-Tu with respect to the 35 remaining inoculation cultures, and characterized using SDS-PAGE gels and eGFP expression. **(b)** Concentration of EF-Tu in OnePot PURE reactions derived from SDS-PAGE gel analysis, as a function of relative volume ratios of the EF-Tu inoculation culture in a co-culture. Each data point represents four biological replicates (mean \pm s.d.). **(c)** *In vitro* eGFP expression activity after 3h plotted against concentration of EF-Tu in OnePot PURE reactions. Measurements on the x-axis represent biological replicates, and y-axis measurements represent four biological replicates with three technical replicates. Error bars represent s.d. **(d)** Coomassie blue-stained SDS-PAGE gel of four OnePot PURE (EF-Tu 47%) replicates. In the right panel, intensities of the different replicates are plotted with molecular weight standards (kDa). **(e)** *In vitro* eGFP expression activity after 3 h plotted against relative inoculation volume ratios of EF-Tu. Each data point represents a single biological replicate with three technical replicates; error bars represent s.d. of the technical replicates. **(f)** Time course of *in vitro* eGFP expression with OnePot PURE (EF-Tu 47%). Each line represents a technical replicate.

the protein components led to only negligible changes in protein expression yield of the PURE system (data not shown). Additionally, over 50% of the HomeMade PURE protein components consists of EF-Tu (Supplementary Table S5). Hence, we decided to optimize our OnePot PURE system with a particular focus on this translation factor.

We varied the relative volume of the EF-Tu inoculating culture with respect to the 35 remaining inoculation cultures to generate ratios of 3%, 17%, 38%, and 47%. The 3% ratio corresponds to 100 μ L of all 36 inoculation cultures, including EF-Tu, being added to the mixed co-culture (Supplementary Table S4). As can be seen from gels and corresponding analysis, larger percentages of the EF-Tu strain in the co-culture led to higher absolute levels of EF-Tu in the OnePot protein system (Fig. 1b, Supplementary Fig. S2, S3). Increased concentrations of EF-Tu also gave rise to higher protein expression yields (Fig. 1c). We could therefore show that it is possible to modify the ratio of an individual PURE protein component simply by varying the initial inoculation ratio of the corresponding strain, and

that the OnePot PURE system gave rise to high protein expression yields.

It has been thought that precise control over the PURE system composition is required to achieve reproducible, and high protein expression yields and it has been suggested that a simple one-pot method would not be a viable option for robustly generating the PURE system [17]. However, we observed that variations in overnight culture densities (Supplementary Fig. S4) did not lead to substantial differences in OnePot PURE protein content (Fig. 1d, Supplementary Fig. S2). We observed high protein expression robustness across four biological replicates, especially for the 38% and 47% EF-Tu formulations, with coefficients of variation (CV) of 8% and 12%, respectively (Fig. 1e, f). In comparison, the CV for a technical replicate of PURExpress and HomeMade PURE were 5% and 12%, respectively. To avoid significant total protein concentration differences across replicates, we adjusted the concentration of the protein mixture to 1.6 mg/mL in the final reaction. This optimal concentration was chosen based on titrations of OnePot PURE (47% EF-Tu) replicate A (Supplementary Fig. S5).

We compared the protein composition of our OnePot PURE system to the commercially available PURExpress (NEB) and our HomeMade PURE system prepared based on the Shimizu protocol with minor adjustments [7]. From gels we determined that the overall composition of the PURExpress and HomeMade PURE systems were quite similar to one another as expected (Fig. 2a). Both PURExpress and HomeMade PURE had a higher relative percentage of EF-Tu and a lower total protein concentration (1 mg/mL for HomeMade PURE) than OnePot PURE. The relative intensities of individual proteins in the OnePot PURE deviated from the PURExpress and HomeMade PURE standards although the protein expression yield of the OnePot PURE system (47% EF-Tu) was similar to PURExpress, 1.6 times higher than our HomeMade PURE and 5 times higher than TraMOS (Fig. 2b).

One of the main factors limiting the use of the PURE system is its high cost. We performed a detailed cost analysis of different PURE systems: two systems prepared from individually purified protein components (PURExpress and HomeMade PURE), as well as two systems prepared from batch cultures and pooled purifications (OnePot and TraMOS) (Fig. 2c, Supplementary Table S1, S2, S3, S6). The commercial PURExpress is the most expensive at a cost of 1.36 USD/ μ L followed by TraMOS (0.96 USD/ μ L), HomeMade PURE (0.36 USD/ μ L), and OnePot PURE (0.09 USD/ μ L). For the HomeMade PURE and TraMOS preparations, cost originates primarily from protein components and ribosomes. The OnePot

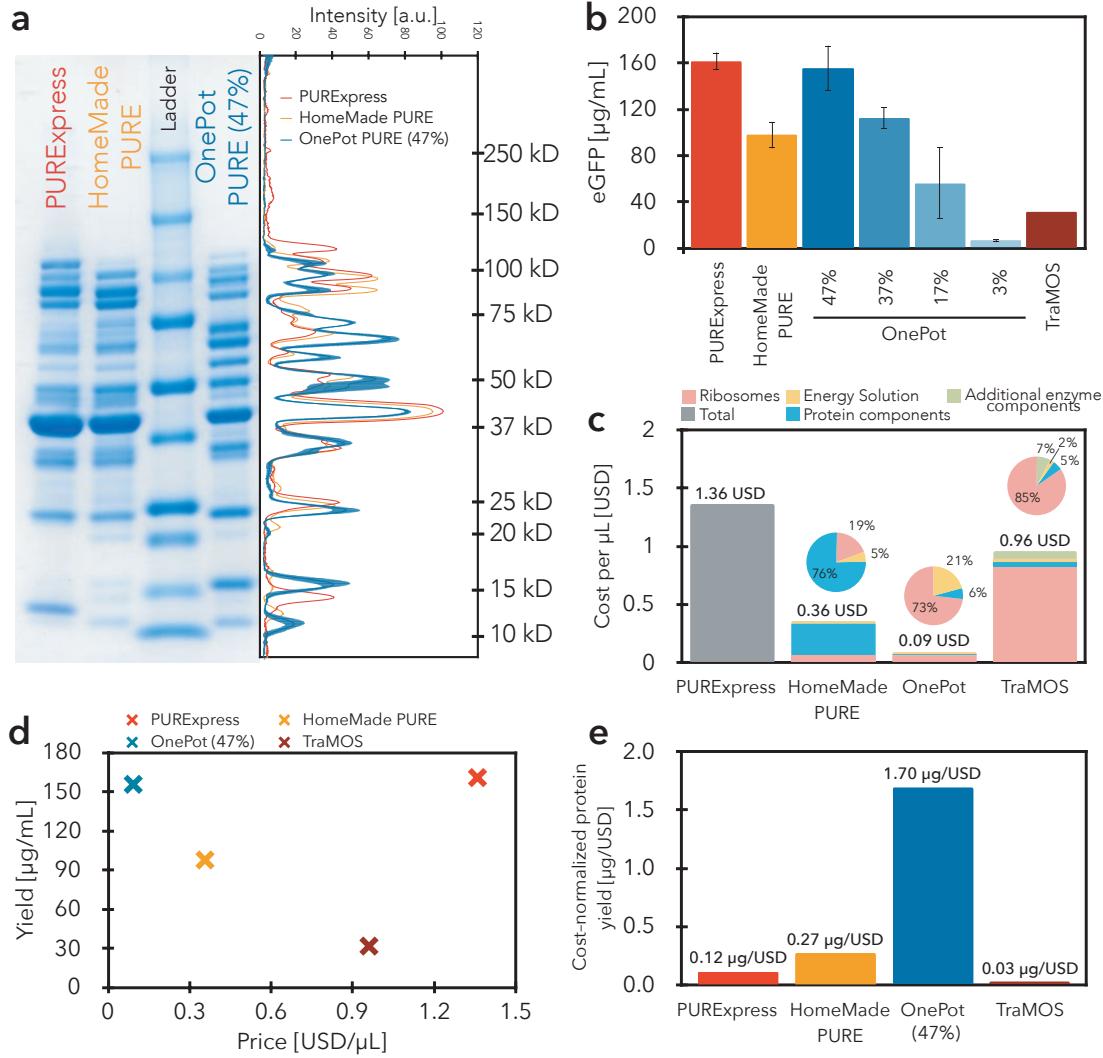


Figure 2: OnePot PURE comparison to existing PURE systems. **(a)** SDS-PAGE gel of PURExpress, HomeMadePURE, OnePot PURE (EF-Tu 47%, replicate A). In the right panel, intensities of different replicates are plotted with molecular weight standards (kDa). **(b)** Comparison of eGFP expression activity (after 3 h) of different PURE systems. The different systems were tested in the same conditions except for TraMOS where the reported value was used [17]. **(c)** Price comparison of the different PURE systems. Calculations are detailed in Supplementary Tables S1, S2, S3 S6. **(d)** Yield of the different PURE systems plotted against their price per μL . Mean values of the eGFP expression yield were plotted. **(e)** Cost-normalized yield of the different PURE systems. The mean value of the eGFP expression yield was used for the calculations.

approach reduces the cost of the protein components to almost negligible levels and relies on ribosome purification to further reduce cost. In-house ribosome purification does not only reduce the price by almost 16-fold as compared to when commercial ribosomes are used, but also allows for higher ribosome concentrations in the PURE system. The standard ribosome purification protocol used in this work is simple and robust. Moreover, in case ultracentrifugation is not accessible, His-tag purification of ribosomes could be a viable alternative [18, 20]. OnePot PURE substantially outperformed all other systems when directly comparing protein synthesis yield and cost per microliter (Fig. 2d) and achieved a cost normalized protein yield of $1.70 \mu\text{g}/\text{USD}$ compared to $0.27 \mu\text{g}/\text{USD}$ for HomeMade PURE, $0.12 \mu\text{g}/\text{USD}$ for PURExpress, and $0.03 \mu\text{g}/\text{USD}$ for TraMOS (Fig. 2e).

We demonstrated that it is possible to robustly produce a highly functional PURE system at low cost using a practical single batch culture and purification approach. The OnePot PURE system described here achieved a protein synthesis yield of $156 \mu\text{g}/\text{mL}$ at a cost of $0.09 \text{ USD}/\mu\text{L}$. At $1.7 \mu\text{g}/\text{USD}$ the cost normalized protein synthesis yield is over a magnitude higher than the commercial PURE system and substantially higher than TraMOS. We also showed that it is possible to adjust and optimize the OnePot PURE system by varying the inoculation fraction of an individual strain. This simple, low-cost, and robust protocol for producing the PURE system should broaden access to the technology and enable new applications which hitherto were not feasible due to the high cost and complexity of producing the PURE system.

MATERIALS AND METHODS

Escherichia coli strains and plasmids

E. coli BL21(DE3) and M15 strains were used for protein expression. All plasmids encoding PURE proteins used in this work were originally obtained from Y. Shimizu (RIKEN Quantitative Biology Center, Japan). Genes coding for MK and PPiase were originally cloned in pET29b vectors with kanamycin resistance. To establish a OnePot system, we used CPEC assembly (Circular Polymerase Extension Cloning) [21] to clone a DNA fragment amplified from pET29b vectors containing MK and PPiase genes as well as the T7 promoter, RBS, and T7 terminator, into a pET21a vector containing ampicillin resistance. The primer

sequences used are listed in Supplementary Table S7. A list of the PURE proteins with their corresponding gene, vector and reference number are given in Supplementary Table S8. *E. coli* A19 (Coli Genetic Stock Center, CGSC#: 5997) was used for ribosome purification. The linear template DNA for *in vitro* eGFP synthesis was initially prepared by extension PCR from a pKT127 plasmid as described [2] and cloned into a pSBlue-1 plasmid. The DNA fragment used for PURE system characterization was amplified from this plasmid by PCR (Supplementary Table S7) and purified by DNA Clean and Concentrator-25 (Zymo Research). DNA was eluted in nuclease-free water instead of elution buffer. For all reactions, a linear template from a single purification was used.

Buffers used for protein and ribosome purification

All buffers used in this work are listed in Supplementary Table S9. All buffers were filtered (Flow Bottle Top Filters, 0.45 μ m aPES membrane) and stored at 4°C. 2-mercaptoethanol was added immediately before use.

OnePot protein preparation

Lysogeny broth (LB) used for OnePot protein component preparation was supplemented with 100 μ g/mL ampicillin and all cultures were grown at 37°C, 260 RPM. To allow for fast and easy inoculation, the different strains were stored as a glycerol stock in a single 96 well microplate. All overnight cultures were inoculated by a 96-well replicator (VP 408FS2AS, V & P Scientific), except for the EF-Tu strain, and grown in 0.3 mL of LB in a deep-well microplate (96 wells, void volume 1.5 mL). The strain expressing EF-Tu was grown in 3 mL of LB in a standard 14 mL culture tube. Overnight cultures (in total 3.6 mL) were used to inoculate 500 mL of LB media in a 1 L baffled flask. The exact composition of the inoculation cultures for different OnePot systems are given in Supplementary Table S4. Cells were grown 2 h before induction with 0.1 mM IPTG for 3 h, then harvested by centrifugation (4,000 RPM, 10 min, 4°C) and stored at -80°C overnight. Cells were resuspended in 7.5 mL buffer A and lysed by sonication on ice (Vibra cell 75186 and probe tip diameter: 6 mm, 4 \times 20s:20s pulse, 70% amplitude). Cell debris was removed by centrifugation (15,000 RPM, 20 min, 4°C). The supernatant was mixed with 2 mL of equilibrated resin, prepared

as described below, and incubated for 3 h, at 4°C. After the incubation, unbound lysate was allowed to flow through the column. The column was washed with 25 mL of a wash buffer (95% buffer A, 5% buffer B) and eluted with 5 mL of elution buffer (10% buffer A, 90% buffer B). Instead of dialysis, buffer exchange was done using a 15 mL Amicon Ultra filter unit with a 3 kDa molecular weight cutoff (Merck). All centrifugation steps were performed at 4,000 RPM and 4°C. The elution fraction was diluted with 25 mL of HT buffer and concentrated to 1 mL (2 × 60 min). The concentrated sample was then diluted with 10 mL of HT buffer, concentrated to 1.5 mL (60 min), and mixed with 1.5 mL of stock buffer B. The protein solution was then concentrated (14,000 RPM, 30 min, 4°C) using a 0.5 mL Amicon Ultra filter unit with a 3 kDa molecular weight cutoff (Merck) and stored at –80°C. Total protein concentration in the OnePot protein mixture was determined using a microplate Bradford protein assay with bovine gamma-globulin as a standard (Bio-Rad). Samples were diluted 1:25 and 5 μL of the diluted sample was mixed with 250 μL of Bradford reagent. Absorbance at 595 nm was measured using a SynergyMX platereader (BioTek). The OnePot protein mixture was then adjusted to a concentration of 12.25 mg/mL.

HomeMade PURE protein preparation

Proteins were prepared by Ni-NTA gravity-flow chromatography. The LB medium used was supplemented with 100 μg/mL of ampicillin and/or 50 μg/mL of kanamycin (Supplementary Table S5), and all cultures were grown at 37°C, 250 RPM. Overnight cultures were grown in 3 mL of LB. Each strain was then individually inoculated in a flask with 2 L of LB. Cells were grown 2 h before induction with 0.1 mM of IPTG for 3 h, then harvested by centrifugation and stored at –80°C overnight. The cells were resuspended in 30 mL of buffer A and lysed by sonication on ice (Vibra cell 75186 and probe tip diameter: 6 mm, 8 × 20s:20s pulse, 70% amplitude). Cell debris was removed by centrifugation (25,000 RCF, 20 min, 4°C). The supernatant was mixed with 2-3 mL of equilibrated resin (described below), and incubated for 1-2 h, at 4°C. After the incubation, unbound lysate was allowed to flow through the column. The column was washed with 30 mL of a wash buffer (95% buffer A, 5% buffer B) and eluted with 15 mL of an elution buffer (10% buffer A, 90% buffer B). The elution fraction was dialysed against HT buffer (2×) and stock buffer and stored at –80°C. Protein concentrations were estimated by absorbance at 280 nm and calculated protein ex-

tinction coefficients. When a higher protein concentration was required, the protein solution was concentrated using a 0.5 mL Amicon Ultra filter unit (Merck).

Ni-NTA resin preparation and regeneration

2 mL IMAC Sepharose 6 FF (GE Healthcare) was pipetted into Econo-Pac chromatography columns (Bio-Rad), and charged with 15 mL of 100 mM nickel sulfate solution. The charged column was washed with 50 mL of DEMI water and equilibrated with 35 mL of buffer A. After protein purification, columns were regenerated with 10 mL of buffer containing 0.2 M EDTA and 0.5 M NaCl, and washed with 30 mL of 0.5 M NaCl, followed by 30 mL of demineralized water, and stored in 20% ethanol at 4°C.

OD600 measurement

OD600 measurements of over-night cultures were measured on a 96-well plate with tenfold dilutions (20 μ L of over-night culture in 180 μ L of LB) using a SynergyMX platereader (BioTek). The background (OD600 of 200 μ L of LB) was subtracted from all samples.

Ribosome purification

Ribosomes were prepared from *E. coli* A19 by hydrophobic interaction chromatography (HIC) and sucrose cushion buffer ultracentrifugation as described previously with slight modifications [22, 23]. *E. coli* A19 strain was grown overnight in 100 mL of LB media at 37°C. 2 \times 30 mL of the overnight cultures was used to inoculate 2 \times 2 L of LB. Cells were grown at 37°C, 250 RPM to exponential phase (3-4 h, OD600 = 0.6-0.8), harvested by centrifugation (4,000 RCF, 20 min, at 4°C), resuspended in 50 mL suspension buffer and stored at -80°C. The resuspended cells were lysed by sonication on ice (Vibra cell 75186 and probe tip diameter: 6 mm, 12 \times 20s:20s pulse, 70% amplitude). The cell debris was removed by centrifugation (20,000 RCF, 20 min, at 4°C). The recovered fraction was mixed with the same amount of high salt suspension buffer. The precipitate was removed by centrifugation (20,000 RCF, 20 min, at 4°C) and the supernatant was filtrated with a GD/X syringe filter membrane (0.45 mm, PVDF, Whatman).

Ribosomes were purified using a 15 mL (3 × 5 mL HiTrap Butyl HP column (GE Healthcare) on Akta Purifier FPLC (GE Healthcare) at 4°C. After the column was equilibrated with 60 mL of buffer C, the prepared lysate solution was loaded onto the column and washed with 45 mL of wash buffer 1 (100% buffer C) followed by 75 mL of wash buffer 2 (80% buffer C, 20% buffer D). Ribosomes were eluted with 60 mL of ribosome elution buffer (50% buffer C, 50% buffer D) followed by 60 mL of final elution buffer (100% buffer D) at a flow rate of 4 mL per minute. All fractions containing ribosomes (absorbance peak at 280 nm during elution with ribosome elution buffer) were pooled together (around 55 mL). The column was recovered by washing with NaOH (1 M) and acetic acid (0.1 M), and stored in 20% ethanol.

14 mL of recovered fraction was overlaid onto 15 mL of cushion buffer in four polycarbonate tubes (void volume: 32 mL). The ribosomes were pelleted by ultracentrifugation (Beckman type SW 32 Ti rotor, 100,000 RCF, 16 h, 4°C). Each transparent ribosome pellet was washed two times with 0.5 mL ribosome buffer and resuspended with a magnetic stirrer in 100 μ M of ribosome buffer. To ensure that all the ribosomes are recovered every tube was washed with 100 μ M ribosome buffer. The recovered solution was concentrated using a 0.5 mL Amicon Ultra filter unit with a 3 kDa molecular weight cutoff (Merck) by centrifugation (14,000 RCF, 10 min, at 4°C). Ribosome concentrations were determined by measuring absorbance at 260 nm of a 1:100 dilution. An absorbance of 10 for the diluted solution corresponded to a 23 μ M concentration of undiluted ribosome solution. Final ribosome solution used for *in vitro* protein synthesis was prepared by diluting the sample to 10 μ M. The usual yield is above 0.75 mL of 10 μ M ribosome solution.

SDS-PAGE gels

Proteins were separated by SDS-PAGE using 15-well 4-20% Mini-PROTEAN TGX Pre-cast Protein Gels (Bio-Rad). Gels were stained using Bio-Safe Coomassie stain (Bio-Rad), scanned using an EPSON Perfection V10 scanner and analyzed with ImageJ.

Energy solution preparation

Energy solution was prepared as described previously with slight modifications [23]. 2.5× energy solution contained 0.75 mM of each amino acid, 29.5 mM of magnesium acetate,

250 mM of potassium glutamate, 5 mM of ATP and GTP, 2.5 mM CTP, UTP, and DTT (Dithiothreitol), 130 U_{A260}/mL of tRNA, 50 mM of creatine phosphate, 0.05 mM of folinic acid, 5 mM of spermidine, and 125 mM of HEPES.

***In vitro* protein expression**

HomeMade or OnePot PURE reactions (5 μ L) were established by mixing 2 μ L of 2.5 \times energy solution, 1.45 μ L linear eGFP DNA template (final concentration: 5 nM), 0.9 μ L of 10 μ M ribosomes (final concentration: 1.8 μ M) and 0.65 μ L of PURE proteins (HomeMade or OnePot solution). PURExpress reactions (5 μ L) were established by mixing 2 μ L of solution A, 1.5 μ L eGFP DNA (final concentration: 5 nM) and 1.5 μ L of solution B. Reactions were incubated at 37°C at constant shaking for 3 h, and measured (excitation: 488 nm, emission: 507 nm) on a SynergyMX platereader (BioTek). Absolute eGFP concentrations were determined from a standard curve (Supplementary Fig. S7).

Cost calculations

To estimate the cost of PURE systems, we analyzed in detail the costs of the different subsets: protein components, ribosomes, and energy solution. The calculation for protein subset costs varies with the type of the system. For the TrAMOS system, the reported cost of 0.052 USD/ μ L was used [17]. For our OnePot system, the cost was estimated based on the calculations given in Supplementary Table S1, with the assumptions that some of the materials can be reused and that four purifications can be done simultaneously in one working day. In the case of the HomeMade PURE system, our estimate was based on the price charged by the EPFL protein expression core facility: 300 USD per 2 L expression culture, which corresponds to our calculation for OnePot PURE of 83 USD per 0.5 L culture (332 USD for 2 L, Supplementary Table S1). Although the total price of this PURE system is high, the total amount of proteins purified is higher as well which can generate at least 40 mL of PURE HomeMade system (based on the volume of the protein limiting the preparation, in our case EF-Tu). Therefore, the price per μ L of HomeMade protein components is 0.27 USD.

Two different possibilities were taken into account in the case of the ribosome subset. In

the first system, commercial ribosomes (Supplementary Table S6) were used for the PURE reactions (TraMOS). In the second system, purified ribosomes were used (HomeMade and OnePot PURE). The cost calculations for purified ribosomes are given in Supplementary Table S2, with the assumptions that some of the materials can be reused and that hands-on time for one purification is a single working day.

The cost calculation for the OnePot energy solution is described in Supplementary Table S3, with the assumption that half a day is necessary for the preparation of 20 mL of energy solution. For the TraMOS energy solution and the additional protein components, the costs were recalculated based on the component's price that would apply for the preparation of the given solutions (Supplementary Table S6). For some of the additional protein components, we were not able to determine the exact protein which was purchased and its amount used, mostly due to a difference in the type of units reported in the paper as compared to the units specified by the supplier. However, we arrived at a very similar cost estimate as given in the original calculation. Furthermore, we assumed that the work required for the solution preparation is taken into account in the purification cost calculation, so we did not consider it.

In the case of PURExpress, the total cost was based on the commercial price. The values used in the cost calculation were derived from experience with the actual experiments while preparing the different subsets. All costs for the different components were based on the prices given in our internal EPFL system when performing the calculation; no delivery costs were taken into account.

ACKNOWLEDGEMENTS

We thank Yoshihiro Shimizu (RIKEN) for kindly supplying plasmids for the PURE system, members of the EPFL Protein Crystallography and Protein Expression Core Facilities, especially Jean Philippe Gaudry for providing us with purified PURE protein components and help with ribosome purification and Nadanai Laohakunakorn for helpful discussions. We also thank Corinna Tuckey and John DeMartino from NEB for providing us with custom PURExpress kits. This work was supported by the École Polytechnique Fédérale de Lausanne.

AUTHOR CONTRIBUTIONS

B.L. performed experiments. B.L. and S.J.M. designed experiments, analyzed data, and wrote the manuscript.

DECLARATION OF INTERESTS

The authors declare no competing interests.

REFERENCES

- [1] Noireaux, V., Bar-Ziv, R. & Libchaber, A. Principles of cell-free genetic circuit assembly. *Proceedings of the National Academy of Sciences of the United States of America* **100**, 12672–12677 (2003).
- [2] Niederholtmeyer, H., Stepanova, V. & Maerkl, S. J. Implementation of cell-free biological networks at steady state. *Proceedings of the National Academy of Sciences of the United States of America* **110**, 15985–15990 (2013).
- [3] Niederholtmeyer, H. *et al.* Rapid cell-free forward engineering of novel genetic ring oscillators. *eLife* **4**, e09771 (2015).
- [4] Karzbrun, E., Tayar, A. M., Noireaux, V. & Bar-Ziv, R. H. Programmable on-chip DNA compartments as artificial cells. *Science (New York, NY)* **345**, 829–832 (2014).
- [5] Goering, A. W. *et al.* In vitro reconstruction of nonribosomal peptide biosynthesis directly from DNA using cell-free protein synthesis. *ACS Synthetic Biology* **6**, 39–44 (2017).
- [6] Nguyen, P. H. B., Wu, Y., Guo, S. & Murray, R. M. Design Space Exploration of the Violacein Pathway in Escherichia coli Based Transcription Translation Cell-Free System (TX-TL). *bioRxiv* 027656 (2016).
- [7] Shimizu, Y. *et al.* Cell-free translation reconstituted with purified components. *Nature Biotechnology* **19**, 751–755 (2001).
- [8] van Nies, P. *et al.* Self-replication of DNA by its encoded proteins in liposome-based synthetic cells. *Nature communications* 1–12 (2018).
- [9] Pardee, K. *et al.* Rapid, Low-Cost Detection of Zika Virus Using Programmable Biomolecular Components. *Cell* **165**, 1255–1266 (2016).

- [10] Pardee, K. *et al.* Portable, On-Demand Biomolecular Manufacturing. *Cell* **167**, 248–254.e12 (2016).
- [11] Huang, A. *et al.* BiobitsTM explorer: A modular synthetic biology education kit. *Science Advances* **4**, eaat5105 (2018).
- [12] Forster, A. C. & Church, G. M. Towards synthesis of a minimal cell. *Molecular systems biology* **2**, 45 (2006).
- [13] Schwille, P. *et al.* MaxSynBio - Avenues towards creating cells from the bottom up. *Angewandte Chemie-International Edition* (2018).
- [14] Li, J., Gu, L., Aach, J. & Church, G. M. Improved Cell-Free RNA and Protein Synthesis System. *PLoS ONE* **9**, e106232–11 (2014).
- [15] Zhang, Y., Huang, Q., Deng, Z., Xu, Y. & Liu, T. Enhancing the efficiency of cell-free protein synthesis system by systematic titration of transcription and translation components. *Biochemical Engineering Journal* **138**, 1–26 (2018).
- [16] Maddalena, L. L. d. *et al.* GreA and GreB Enhance Expression of Escherichia coliRNA Polymerase Promoters in a Reconstituted Transcription–Translation System. *ACS Synthetic Biology* **5**, 929–935 (2016).
- [17] Villarreal, F. *et al.* Synthetic microbial consortia enable rapid assembly of pure translation machinery. *Nature Chemical Biology* **14**, 29–35 (2017).
- [18] Wang, H. H. *et al.* Multiplexed in vivo His-tagging of enzyme pathways for in vitro single-pot multienzyme catalysis. *ACS Synthetic Biology* **1**, 43–52 (2012).
- [19] Li, G.-W., Burkhardt, D., Gross, C. & Weissman, J. S. Quantifying absolute protein synthesis rates reveals principles underlying allocation of cellular resources. *Cell* **157**, 624–635 (2014).
- [20] Ederth, J., Mandava, C. S., Dasgupta, S. & Sanyal, S. A single-step method for purification of active His-tagged ribosomes from a genetically engineered Escherichia coli. *Nucleic acids research* **37**, e15 (2009).
- [21] Quan, J. & Tian, J. Circular polymerase extension cloning of complex gene libraries and pathways. *PLoS ONE* **4**, e6441 (2009).
- [22] Horiya, S., Bailey, J. K. & Krauss, I. J. Directed Evolution of Glycopeptides Using mRNA Display. In *Methods in Enzymology*, 83–141 (Elsevier, 2017).
- [23] Shimizu, Y. & Ueda, T. PURE technology. *Methods Mol Biol* **607**, 11–21 (2010).

SUPPLEMENTAL FIGURES AND TABLES

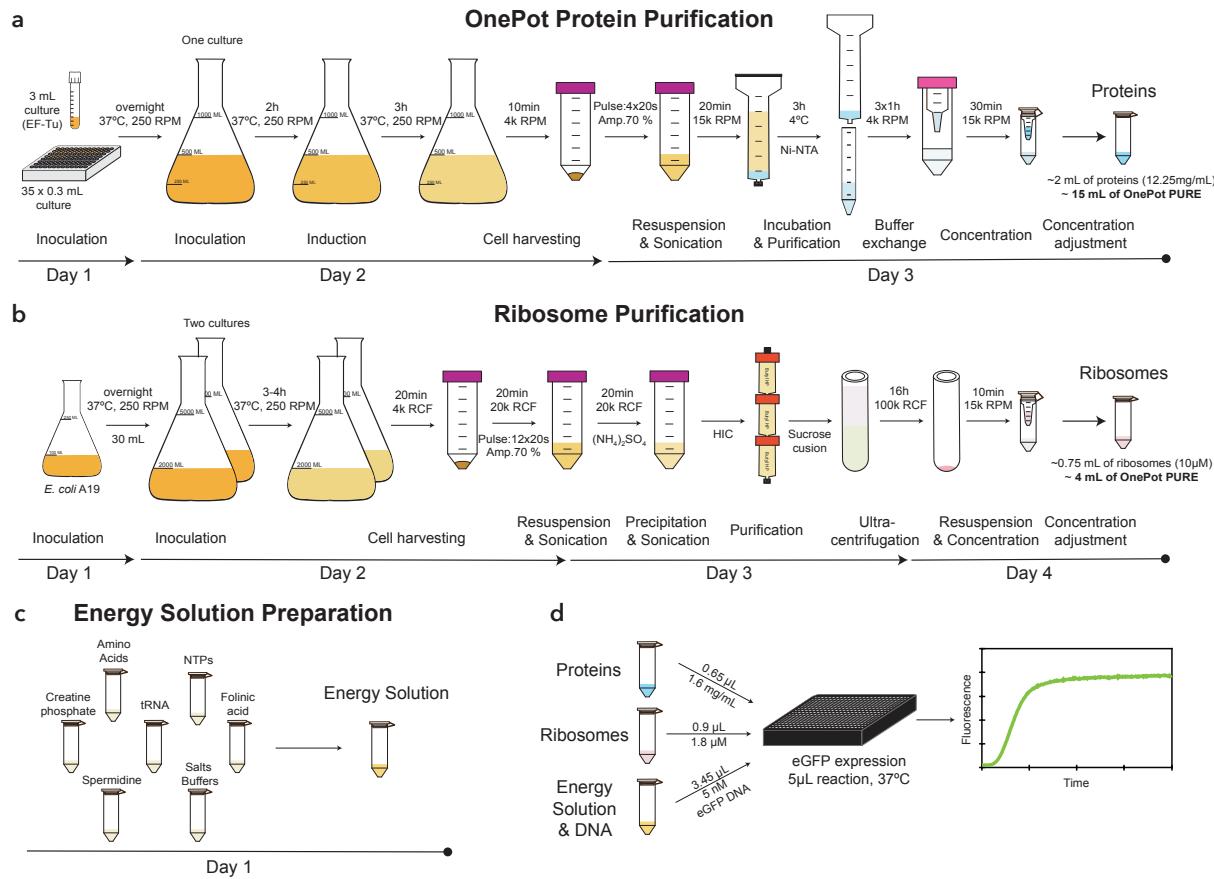


Figure S1: Schematics depicting all steps of the OnePot PURE production. **(a)** Protein purification, **(b)** ribosome purification, and **(c)** energy solution preparation steps. The description of the different steps as well as the day on which they are performed are indicated below the schematics. **(d)** Composition of the OnePot PURE reaction. Two numbers are given for each subset, the volume required for a 5 μ L reaction and the component concentration in the reaction.

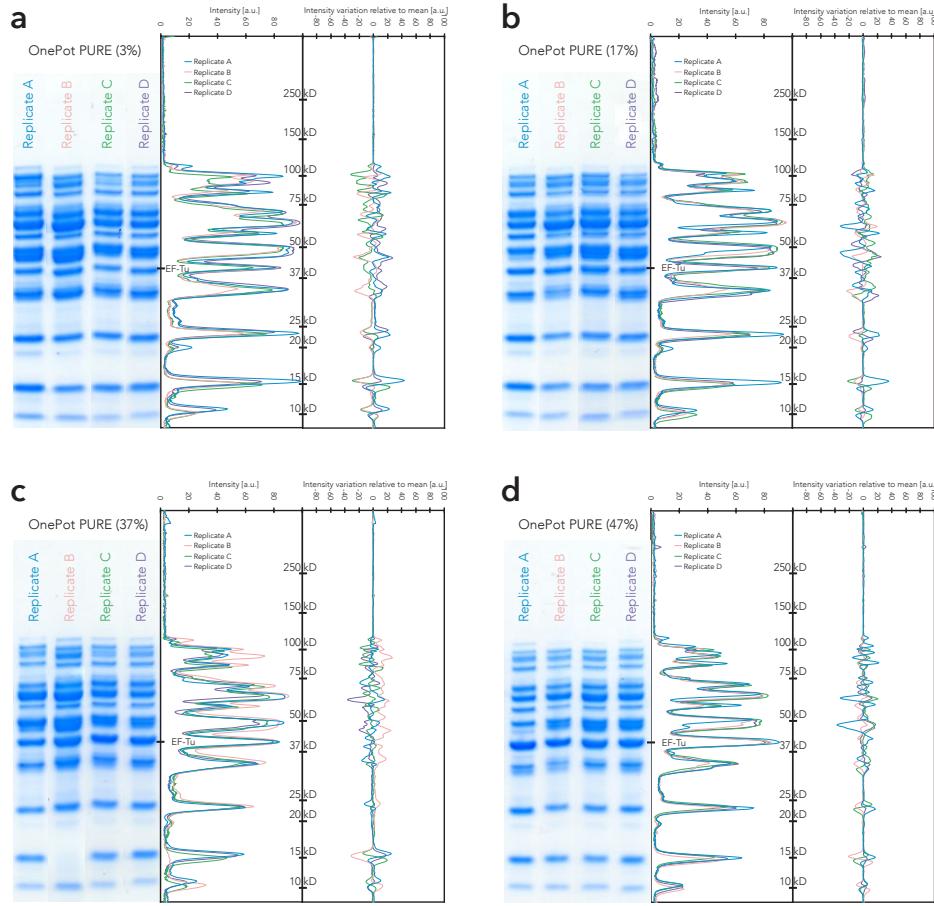


Figure S2: Coomassie blue stained SDS-PAGE gels of the four OnePot PURE formulations. **(a)** 3% EF-Tu, **(b)** 17% EF-Tu, **(c)** 37% EF-Tu, and **(d)** 47% EF-Tu. In the panels to the left of the gels, intensities of the different replicates are plotted with molecular weight standards (kDa). On the right the intensity variations relative to the inter-replicate mean is shown.

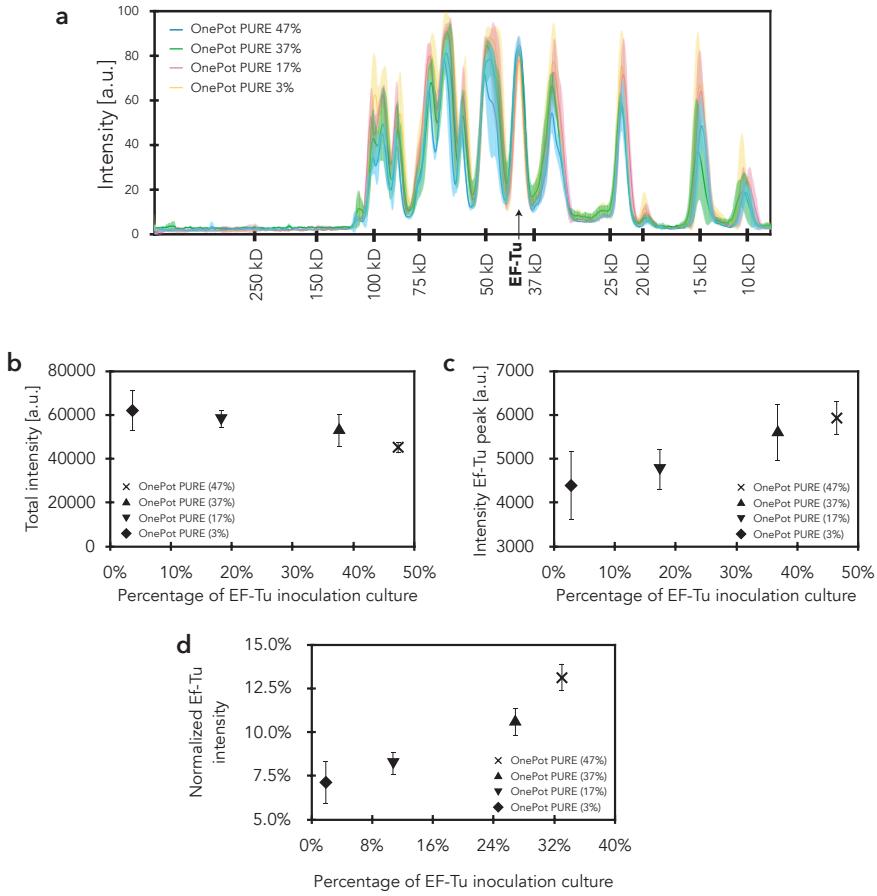


Figure S3: EF-Tu analysis. **(a)** Mean intensities of the different OnePot systems are plotted against molecular weight standards (kDa); the shaded regions represent the s.d. of the four biological replicates. **(b)** Total intensity of all protein bands as a function of EF-Tu clone inoculation percentage. **(c)** The integrated intensity of the EF-Tu peak from SDS-PAGE gel analysis as a function of EF-Tu clone inoculation percentage. **(d)** The normalised EF-Tu intensity (integrated EF-Tu peak intensity / total protein intensity) as a function of EF-Tu clone inoculation percentage. **(b) - (d)** Each data point represents four biological replicates (mean \pm s.d.)

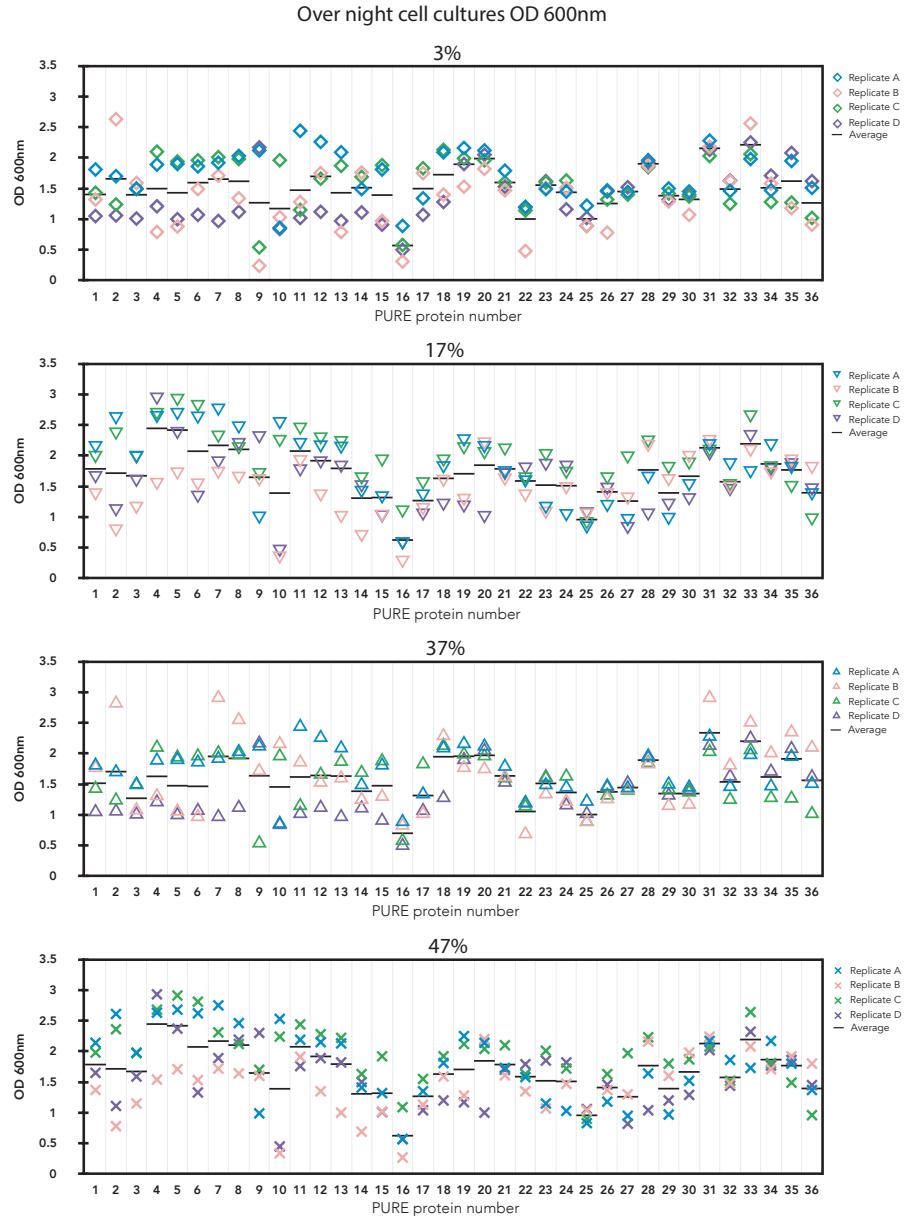


Figure S4: OD₆₀₀ measurement of overnight cultures used for inoculation of the mixed culture.

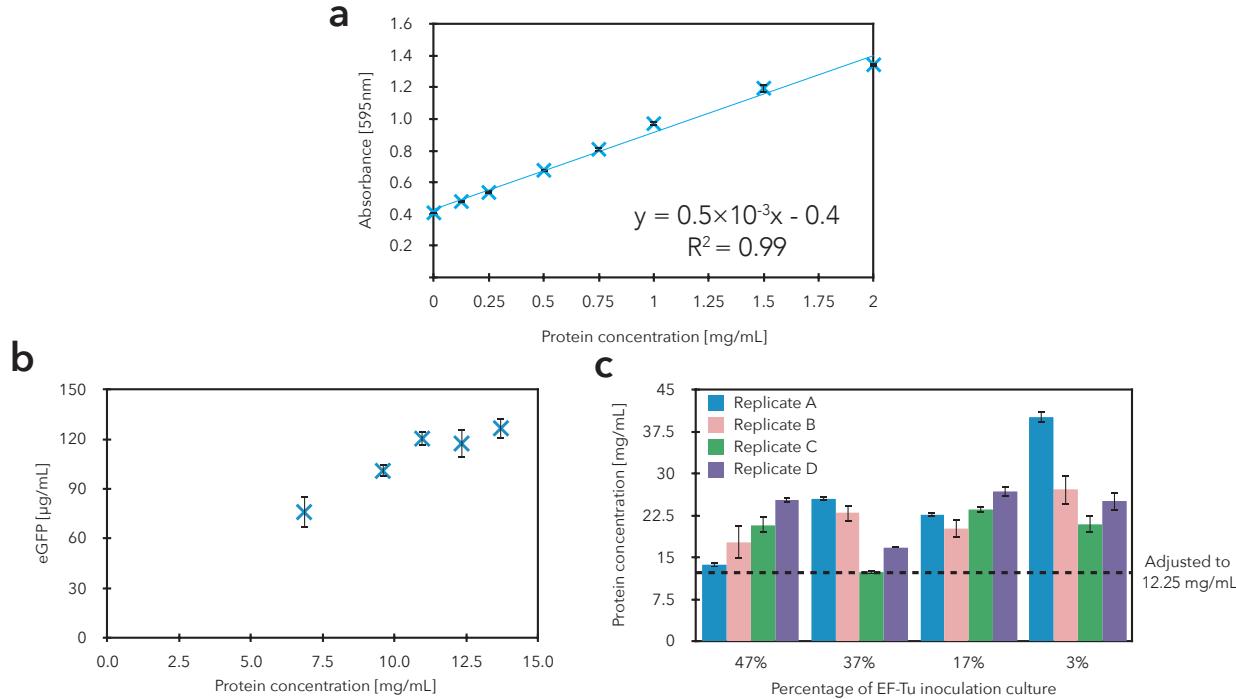


Figure S5: Protein concentration calibrations and adjustments. **(a)** Bradford assay standard calibration curve for protein concentration. The standard curve was produced by measuring the absorbance at 595 nm of prediluted bovine γ -globulin standards. Data are shown as mean \pm s.d. ($n = 3$). Linear fit errors were not propagated as they were negligible compared to experimental errors. **(b)** eGFP expression as a function of protein concentrations in the protein subset of OnePot PURE (47%) replicate A (7.7 \times concentration in the final reaction). Each point represents at least two replicates; data are shown as mean \pm s.d. **(c)** The concentrations of all OnePot protein subsets and their replicates after purification. Each bar represents two independent measurements in technical duplicate. Data are shown as mean \pm s.d. The dotted line represents concentration (12.25 mg/mL, which is equal to 1.6 mg/mL in the final PURE reaction) to which all reactions were adjusted to.

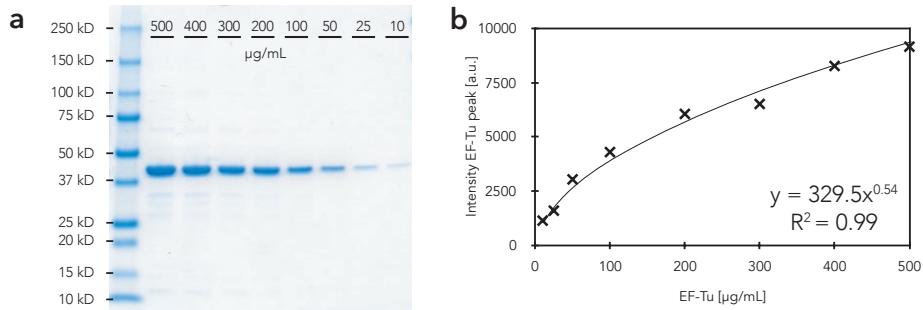


Figure S6: Standard calibration curve for EF-Tu protein concentration. **(a)** Coomassie blue stained SDS-PAGE gels of different concentrations of EF-Tu. **(b)** The standard curve was produced by measuring the integrated intensity of the EF-Tu peak at different EF-Tu concentrations. The reference EF-Tu concentration was determined by absorbance measurement at 280 nm.

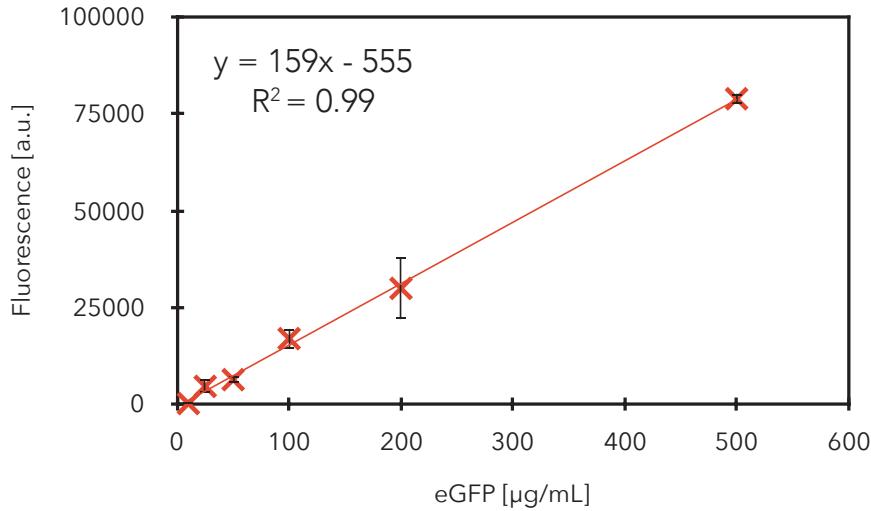


Figure S7: Standard calibration curve for eGFP. The standard curve was produced by measuring the fluorescence over 60 min for different eGFP (TP790050, AMS Biotechnology) concentrations in PBS on a plate reader with the same settings as for *in vitro* expression. Excitation and emission wavelengths were 488 nm and 507 nm, respectively. Experiments were performed in triplicates. Fluorescence measurements for the first 20 min were not considered. Data are shown as mean \pm s.d. ($n = 3$). Linear fit errors were not propagated as they were negligible compared to experimental errors.

TABLE S1: OnePot protocol cost and time estimate

OnePot Protein Purification

Description	Catalog Number	Company	Amount	Price [USD]	Amount per purification	Price per purification [USD]	Note
LB media	A0954	PanReac AppliChem	2,500 g	163	15 g	0.98	
IMAC Sepharose® 6 Fast Flow	17-0921-07	GE Healthcare	25 mL	208	0.5 mL	4.17	2mL per purification (reused at least for 4 purifications)
Econo-Pac Chromatography Columns	7321010	Bio-Rad Laboratories	50 pcs	382	0.25 pcs	1.91	1 per purification (reused at least for 4 purifications)
Nickel Sulfate	15414469	Alfa Aesar	100 mL	47	3 mL	1.41	
Buffers						4.30	
AMICON ULTRA 15ML - 3 KDa	UFC900324	Merck Millipore	24 pcs	248	1 pcs	10.33	
AMICON ULTRA 0.5ML - 3 KDa	Merck Millipore	Merck Millipore	24 pcs	112	1 pcs	4.69	
Additional Lab supplies (pipets, tubes)						5.00	
Work			1	200	0.25	50.00	4 purifications can be done at the same time
Total price per one purification						82.8	
Amount of PURE from single purification			15 mL	Price per 1µL		0.006	

Protein Purification	Total time	Active time
Day 1	14h	2h 30m
Inoculation + Cell Growth	12h	30m
Preparation of media and buffers	2h	2h
Day 2	6h 30m	1h 30m
Inoculation + Cell Growth	2h 20m	20m
Induction + Cell Growth	3h 10m	10m
Centrifugation	1h	1h
Day 3	9h	4h
Preparation of columns	40m	40m
Cell sonication and centrifugation	40m	40m
Purification	3h 30m	1h
Buffer exchange	3h 30m	1h
Concentration	40m	40m
Total	29h 30m	8h

TABLE S2: Ribosome protocol cost and time estimate

Ribosomes Purification

Description	Catalog Number	Company	Amount	Price	Amount per purification	Price per purification [USD]	Note
LB media	A0954	PanReac AppliChem	2500 g	CHF163.27	100 g	6.53	
HiTrap Butyl HP Column	28411005	GE Healthcare	5 pcs	CHF322.73	0.2 pcs	12.91	2-3 per purification (reused for multiple purifications)
Thickwall Polycarbonate Tube	355631	Beckman	25 pcs	CHF274.63	1 pcs	10.99	reused for multiple purifications
Whatman® GD/X syringe filters	WHA68722504	GE Whatman	50 pcs	CHF303.68	1 pcs	6.07	
Buffers						22.63	
AMICON ULTRA 0.5ML - 3 KDa	Merck Millipore	Merck Millipore	24 pcs	CHF112.45	1 pcs	4.69	
Additional Lab supplies (pipets, tubes)						5.00	
Work			1	CHF200.00	1	200.00	
Total price per one purification						268.8	
Amount of PURE from single purification				4 mL	Price per 1µL		0.07

Ribosome Purification	Total time	Active time
Day 1	14h	2h 10m
Inoculation + Cell Growth	12h	10m
Preparation of media and buffers	2h	2h
Day 2	5h	1h 10m
Inoculation + Cell Growth	4h	10m
Centrifugation	1h	1h
Day 3	21h	3h 40m
Preparation and cleaning of columns	2h	1h
Cell sonication and centrifugation	1h 30m	40m
Purification	1h 30m	1h 30m
Ultracentrifugation	16h	30m
Day 3	1h	1h
Resuspension, Concentration	1h	1h
Total	41h	8h

TABLE S3: Energy solution cost estimates

OnePot Energy Solution

Compound	Catalog Number	Company	Amount	Price [USD]	Amount per 1uL reaction	Price per 1uL reaction [USD]
Amino acids	LAA21-1KT	Sigma-Aldrich	1 g	464	0.04 µg	0.000020
Magnesium acetate	M0631	Sigma-Aldrich	100 g	20	2.5 µg	0.000001
Potassium glutamate	49601	Sigma-Aldrich	500 g	120	20 µg	0.000005
DTT	sc-29089B	SantaCruz Biotech	10 g	139	0.15 µg	0.000002
ATP	R0481	ThermoFisher	250 uL	36	20 nL	0.002844
GTP	R0481	ThermoFisher	250 uL	36	20 nL	0.002844
CTP	R0481	ThermoFisher	250 uL	36	10 nL	0.001422
UTP	R0481	ThermoFisher	250 uL	36	10 nL	0.001422
tRNA	10109541001	Roche	500 uL	244	17 nL	0.008186
Creatine phosphate	27920	Sigma-Aldrich	1 g	38	7 µg	0.000248
Folinic acid	F7876	Sigma-Aldrich	25 g	84	0.01 µg	0.000000
Spermidine	S2626	Sigma-Aldrich	1 g	34	0.29 µg	0.000010
HEPES	H0887-100ML	Sigma-Aldrich	100 mL	39	50 nL	0.000020
Work			20 mL	100	0.4 uL	0.002000
Total Price per 1µL						0.019

TABLE S4: OnePot inoculation culture volumes

Number	Protein	Vector	Strain	OnePot (3%)		OnePot (17%)		OnePot (37%)		OnePot (47%)	
Amount of inoculation culture				µL	%	µL	%	µL	%	µL	%
1	AlaRS	pQE30	M15	100	2.8%	85	2.4%	65	1.8%	55	1.5%
2	ArgRS	pET16b	BL21(DE3)	100	2.8%	85	2.4%	65	1.8%	55	1.5%
3	AsnRS	pQE30	M15	100	2.8%	85	2.4%	65	1.8%	55	1.5%
4	AspRS	pET21a	BL21(DE3)	100	2.8%	85	2.4%	65	1.8%	55	1.5%
5	CysRS	pET21a	BL21(DE3)	100	2.8%	85	2.4%	65	1.8%	55	1.5%
6	GlnRS	pET21a	BL21(DE3)	100	2.8%	85	2.4%	65	1.8%	55	1.5%
7	GluRS	pET21a	BL21(DE3)	100	2.8%	85	2.4%	65	1.8%	55	1.5%
8	GlyRS	pET21a	BL21(DE3)	100	2.8%	85	2.4%	65	1.8%	55	1.5%
9	HisRS	pET21a	BL21(DE3)	100	2.8%	85	2.4%	65	1.8%	55	1.5%
10	IleRS	pET21a	BL21(DE3)	100	2.8%	85	2.4%	65	1.8%	55	1.5%
11	LeuRS	pET21a	BL21(DE3)	100	2.8%	85	2.4%	65	1.8%	55	1.5%
12	LysRS	pET21a	BL21(DE3)	100	2.8%	85	2.4%	65	1.8%	55	1.5%
13	MetRS	pET21a	BL21(DE3)	100	2.8%	85	2.4%	65	1.8%	55	1.5%
14	PheRS	pQE30	M15	100	2.8%	85	2.4%	65	1.8%	55	1.5%
15	ProRS	pET21a	BL21(DE3)	100	2.8%	85	2.4%	65	1.8%	55	1.5%
16	SerRS	pET21a	BL21(DE3)	100	2.8%	85	2.4%	65	1.8%	55	1.5%
17	ThrRS	pQE30	M15	100	2.8%	85	2.4%	65	1.8%	55	1.5%
18	TrpRS	pET21a	BL21(DE3)	100	2.8%	85	2.4%	65	1.8%	55	1.5%
19	TyrRS	pET21a	BL21(DE3)	100	2.8%	85	2.4%	65	1.8%	55	1.5%
20	ValRS	pET21a	BL21(DE3)	100	2.8%	85	2.4%	65	1.8%	55	1.5%
21	IF1	pQE30	M15	100	2.8%	85	2.4%	65	1.8%	55	1.5%
22	IF2	pQE30	M15	100	2.8%	85	2.4%	65	1.8%	55	1.5%
23	IF3	pQE30	M15	100	2.8%	85	2.4%	65	1.8%	55	1.5%
24	EF-G	pQE60	M15	100	2.8%	85	2.4%	65	1.8%	55	1.5%
25	EF-Tu	pQE60	M15	100	2.8%	625	17.4%	1325	36.8%	1675	46.5%
26	EF-Ts	pQE60	M15	100	2.8%	85	2.4%	65	1.8%	55	1.5%
27	RF1	pQE30	M15	100	2.8%	85	2.4%	65	1.8%	55	1.5%
28	RF2	pET15b	BL21(DE3)	100	2.8%	85	2.4%	65	1.8%	55	1.5%
29	RF3	pQE30	M15	100	2.8%	85	2.4%	65	1.8%	55	1.5%
30	RRF	pQE60	M15	100	2.8%	85	2.4%	65	1.8%	55	1.5%
31	MTF	pET21a	BL21(DE3)	100	2.8%	85	2.4%	65	1.8%	55	1.5%
32	CK	pQE30	M15	100	2.8%	85	2.4%	65	1.8%	55	1.5%
33	MK	pET21a	BL21(DE3)	100	2.8%	85	2.4%	65	1.8%	55	1.5%
34	NDK	pQE30	M15	100	2.8%	85	2.4%	65	1.8%	55	1.5%
35	PPiase	pET21a	BL21(DE3)	100	2.8%	85	2.4%	65	1.8%	55	1.5%
36	T7 RNAP	pQE30	M15	100	2.8%	85	2.4%	65	1.8%	55	1.5%
Total amount of inoculation culture				3600	100%	3600	100%	3600	100%	3600	100%

TABLE S5: HomeMade PURE protein concentrations

Number	Protein	Vector	Antibiotic	Strain	HomeMade PURE	
					Final Concentration in reaction [µg/mL]	Concentration in PURE protein solution [µg/mL]
1	AlaRS	pQE30	Amp, Kan	M15	70	538
2	ArgRS	pET16b	Amp	BL21(DE3)	2	15
3	AsnRS	pQE30	Amp, Kan	M15	22	169
4	AspRS	pET21a	Amp	BL21(DE3)	8	62
5	CysRS	pET21a	Amp	BL21(DE3)	1	9
6	GlnRS	pET21a	Amp	BL21(DE3)	4	29
7	GluRS	pET21a	Amp	BL21(DE3)	13	97
8	GlyRS	pET21a	Amp	BL21(DE3)	10	74
9	HisRS	pET21a	Amp	BL21(DE3)	1	6
10	IleRS	pET21a	Amp	BL21(DE3)	40	308
11	LeuRS	pET21a	Amp	BL21(DE3)	4	31
12	LysRS	pET21a	Amp	BL21(DE3)	6	49
13	MetRS	pET21a	Amp	BL21(DE3)	2	18
14	PheRS	pQE30	Amp, Kan	M15	17	131
15	ProRS	pET21a	Amp	BL21(DE3)	10	77
16	SerRS	pET21a	Amp	BL21(DE3)	2	15
17	ThrRS	pQE30	Amp, Kan	M15	6	48
18	TrpRS	pET21a	Amp	BL21(DE3)	6	48
19	TyrRS	pET21a	Amp	BL21(DE3)	1	5
20	ValRS	pET21a	Amp	BL21(DE3)	2	14
21	IF1	pQE30	Amp, Kan	M15	10	77
22	IF2	pQE30	Amp, Kan	M15	40	308
23	IF3	pQE30	Amp, Kan	M15	10	77
24	EF-G	pQE60	Amp, Kan	M15	50	385
25	EF-Tu	pQE60	Amp, Kan	M15	500	3846
26	EF-Ts	pQE60	Amp, Kan	M15	50	385
27	RF1	pQE30	Amp, Kan	M15	10	77
28	RF2	pET15b	Amp	BL21(DE3)	10	77
29	RF3	pQE30	Amp, Kan	M15	10	77
30	RRF	pQE60	Amp, Kan	M15	10	77
31	MTF	pET21a	Amp	BL21(DE3)	20	154
32	CK	pQE30	Amp, Kan	M15	4	31
33	MK	pET29b	Kan	BL21(DE3)	3	23
34	NDK	pQE30	Amp, Kan	M15	1	8
35	PPiase	pET29b	Kan	BL21(DE3)	1	8
36	T7 RNAP	pQE30	Amp, Kan	M15	10	77
Total protein concentration [µg/mL]					966	7428

TABLE S6: Energy solution and ribosome cost estimates for TraMOS

TraMOS Energy Solution

Compound	Catalog Number	Company	Amount	Price [USD]	Amount per 1uL reaction	Price per 1uL reaction [USD]
Amino acids	LAA21-1KT	Sigma-Aldrich	1 g	464	0.39 µg	0.000181
Magnesium acetate	M0631	Sigma-Aldrich	100 g	20	3.9 µg	0.000001
Potassium glutamate	49601	Sigma-Aldrich	500 g	120	71 µg	0.000017
DTT	sc-29089B	SantaCruz Biotech	10 g	139	0.77 µg	0.000011
ATP	R0481	ThermoFisher	250 uL	36	38 nL	0.005333
GTP	R0481	ThermoFisher	250 uL	36	25 nL	0.003555
CTP	R0481	ThermoFisher	250 uL	36	13 nL	0.001778
UTP	R0481	ThermoFisher	250 uL	36	13 nL	0.001778
tRNA	10109541001	Roche	500 uL	244	17 nL	0.008501
Creatine phosphate	27920	Sigma-Aldrich	1 g	38	16 µg	0.000620
Folinic acid	F7876	Sigma-Aldrich	25 g	84	0.03 µg	0.000000
Spermidine	S2626	Sigma-Aldrich	1 g	34	0.15 µg	0.000005
HEPES	H0887-100ML	Sigma-Aldrich	100 mL	39	50 nL	0.000020
Total Price per 1µL						0.022

TraMOS Additional Enzymes

Compound	Catalog Number	Company	Amount	Price [USD]	Amount per 1uL reaction	Price per 1uL reaction [USD]
BSA	A3912	Sigma-Aldrich	100 g	229	1 µg	0.000002
Creatine kinase	10127566001	Roche	100 mg	142	81 ng	0.000115
Myokinase	M3003	Sigma-Aldrich	0.3 mg	63	50 ng	0.010425
Diphosphonucleotide kinase	N2635	Sigma-Aldrich	0.1 mg	233	4.1 ng	0.009517
T7 RNAP	M0251S	New England Biolabs	100 uL	96	40 nL	0.038360
RNase inhibitor	M0314S	New England Biolabs	75 uL	101	10 nL	0.013413
Total Price per 1µL						0.072

Ribosomes

Compound	Catalog Number	Company	Amount	Price [USD]	Amount per 1uL reaction	Price per 1uL reaction [USD]
Ribosomes	M0314S	New England Biolabs	1 mg	250	3 µg	0.815163
Total Price per 1µL						0.815

TABLE S7: DNA sequences
Primers used for CPEC

	Forward primers	Reverse primers
Primers used for amplification DNA fragment from pET29b vectors	5'-GCGTCCCATTGCCAATC-3'	5'-GCGTCCCATTGCCAATC-3'
Primers used for amplification DNA fragment from pET21a vectors	5'-CCATTCTTCGGCGG-3'	5'-CTGAAAGGAGGAACATATCCGGATTGG-3'

Reporter DNA sequence and amplification primers

	DNA sequence	Amplification Primers
Reporter eGFP linear DNA fragment	gatcttaaggctagact taatacgact caactatagg ggagaccacaacgggttcccttagaaataattttgtttaacttaag aaggag aaaaaaaaATGCTAAAGGTGAA GAATTATTCACTGGTTGCCCCAATTITGTTGAATTAGGGTGTATGTTAAATGGTCAACAATTITGTTCTCCGGTAAGGTGAA GTGATGCTACTTACGGTAAATTGACCTTAAATTGTTACTGGTAAATTGGCAGTTCCATGCCAACCTTAGTCACACTTTA ACTTATGGTTCAATGTTTCTAGATACCCAGATCATGAACAAACATGACTTTCAAGTCGGCCATGCCAGAAGGTATGTTCA AGAAAGAACTATTTTCAAAGATGACGGTAACTACAAGCAGAGCTGAAGTCAGTTGAAGGTGATACTTAGTTAATAGAAT CGAATTTAAAGGTATTGTTAAAGAAAGATGTAACATTGTTAGTCACAAATTGGAAATACAACATAAACTCTCACATGTTACATCAT GGCTGACAAACAAAAGGATGGTCAAGTTAACTTCAAATTTAGGACAACATGAAAGATGGTTCTGTTCAATTAGCTGACCATTAT CAACAAAATACTCAAATGGTGTGGCCAGTCTGTTACAGACACACCAATTACTTATCCACTCAATCTGCCATTACAAAGATCAA ACGAAAGAGAGACACATGGTCTGTTAGAATTGTTACTGCTGCTGGTATTACCATGGTATGGATGAATTGACAATAAaacgcac taggcgtctacgcctgtactggaaaacaaaaacccaaaaacaaaaactgagccccattggatcgatggaggactctatcaaaaaaaaaaaaaaaaaaaaa aactgataaccccttgggcccttaaacgggtcttgagggggttttttg	5'- GATCTTAAGGCTAGAGTACTA ATACGACTCACTATAGGGAG ACC-3'
Blue	T7 promoter	
Red	RBS	
Green	eGFP	
Bold	T7 terminator	

TABLE S8: PURE protein list

Number	Protein	Protein name	Gene	Organism	Vector used for Home Made PURE	Vector used for OnePot PURE	Expression Strain
1	AlaRS	Alanyl-tRNA synthetase	alaS	E. coli	pQE30	pQE31	M15
2	ArgRS	Arginyl-tRNA synthetase	argS	E. coli	pET16b	pET16b	BL21(DE3)
3	AsnRS	Asparaginyl-tRNA synthetase	asnS	E. coli	pQE30	pQE31	M15
4	AspRS	Aspartate-tRNA synthetase	aspS	E. coli	pET21a	pET21a	BL21(DE3)
5	CysRS	Cysteinyl-tRNA synthetase	cysS	E. coli	pET21a	pET21a	BL21(DE3)
6	GlnRS	Glutaminyl-tRNA synthetase	glnS	E. coli	pET21a	pET21a	BL21(DE3)
7	GluRS	Glutamyl-tRNA synthetase	gltX	E. coli	pET21a	pET21a	BL21(DE3)
8	GlyRS	Glycyl-tRNA synthetase	glyQ & glyS	E. coli	pET21a	pET21a	BL21(DE3)
9	HisRS	Histidyl-tRNA synthetase	hisS	E. coli	pET21a	pET21a	BL21(DE3)
10	IleRS	Isoleucyl-tRNA synthetase	ileS	E. coli	pET21a	pET21a	BL21(DE3)
11	LeuRS	Leucyl-tRNA synthetase	leuS	E. coli	pET21a	pET21a	BL21(DE3)
12	LysRS	Lysyl-tRNA synthetase	lysS	E. coli	pET21a	pET21a	BL21(DE3)
13	MetRS	Methionine-tRNA ligase	metG	E. coli	pET21a	pET21a	BL21(DE3)
14	PheRS	Phenylalanyl-tRNA synthetase	pheT & pheS	E. coli	pQE30	pQE31	M15
15	ProRS	Prolyl-tRNA synthetase	proS	E. coli	pET21a	pET21a	BL21(DE3)
16	SerRS	Seryl-tRNA synthetase	serS	E. coli	pET21a	pET21a	BL21(DE3)
17	ThrRS	Threonyl-tRNA synthetase	thrS	E. coli	pQE30	pQE31	M15
18	TrpRS	Tryptophanyl-tRNA synthetase	trpS	E. coli	pET21a	pET21a	BL21(DE3)
19	TyrRS	Tyrosyl-tRNA synthetase	tyrS	E. coli	pET21a	pET21a	BL21(DE3)
20	ValRS	Valyl-tRNA synthetase	valS	E. coli	pET21a	pET21a	BL21(DE3)
21	IF1	Initiation factor 1	infA	E. coli	pQE30	pQE31	M15
22	IF2	Initiation factor 2	infB	E. coli	pQE30	pQE31	M15
23	IF3	Initiation factor 3	infC	E. coli	pQE30	pQE31	M15
24	EF-G	Elongation factor G	fusA	E. coli	pQE60	pQE61	M15
25	EF-Tu	Elongation factor Tu	tufB	E. coli	pQE60	pQE61	M15
26	EF-Ts	Elongation factor Ts	tsf	E. coli	pQE60	pQE61	M15
27	RF1	Release factor 1	prfA	E. coli	pQE30	pQE31	M15
28	RF2	Release factor 2	prfB	E. coli	pET15b	pET15b	BL21(DE3)
29	RF3	Release factor 4	prfC	E. coli	pQE30	pQE31	M15
30	RRF	Ribosome recycling factor	frr	E. coli	pQE60	pQE61	M15
31	MTF	Methionyl-tRNA formyltransferase	fmt	E. coli	pET21a	pET21a	BL21(DE3)
32	CK	Creatine kinase	CKM	Chicken	pQE30	pQE31	M15
33	MK	Adenylate kinase (Myokinase)	AK1	Chicken	pET29b	pET21a	BL21(DE3)
34	NDK	Nucleotide diphosphate kinase	ndk	E. coli	pQE30	pQE31	M15
35	PPiase	Inorganic pyrophosphatase	IPP1	<i>Saccharomyces cerevisiae</i> <i>Enterobacteria phage T7</i>	pET29b	pET21a	BL21(DE3)
36	T7 RNAP	T7 RNA polymerase	1		pQE30	pQE31	M15

TABLE S9: Buffers and energy solution

Bufferes for protein purification

Compound	Catalog number	Company	Buffer A	Buffer B	Buffer HT	Stock buffer	Stock buffer B	Note
			mM	mM	mM	mM	mM	
HEPES	H0887-100ML	Sigma-Aldrich	50	50	50	50	50	pH = 7.6, KOH
Ammonium chloride	09718-250G	Sigma-Aldrich	1000					
Magnesium chloride	63020-1L	Honeywell Fluka	10	10	10	10	10	
Potassium chloride	P5405-1KG	Sigma-Aldrich		100	100	100	100	
Imidasol	I2399	Sigma-Aldrich		500				pH = 7.6, KOH
Glycerol	G7757-1L	Sigma-Aldrich				30%	60%	
β-mercaptoethanol	M6250-100ML	Sigma-Aldrich	7	7	7	7	7	

Bufferes for ribosomes purification

Compound	Catalog number	Company	Suspension buffer	Suspension buffer high salt	Buffer C	Buffer D	Cusion buffer	Ribosome buffer	Note
			mM	mM	mM	mM	mM	mM	
HEPES	H0887-100ML	Sigma-Aldrich	10	10	20	20	20	20	pH = 7.6, KOH
Magnesium acetate	M0631	Sigma-Aldrich	10	10	10	10	10	6	
Potassium chloride	P5405-1KG	Sigma-Aldrich	50	50				30	
Ammonium chloride	09718-250G	Sigma-Aldrich					30		
Ammonium sulfate	A4418	Sigma-Aldrich		3000	1500				pH = 7.6, KOH
Sucrose	84097	Sigma-Aldrich					30%		
β-mercaptoethanol	M6250-100ML	Sigma-Aldrich	7	7	7	7	7	7	

Energy Solution

Compound	Catalog number	Company	Concentration in reaction	Concentration in subset (2.5x)	Units
Amino acids	LAA21-1KT	Sigma-Aldrich	0.3	0.75	mM
Magnesium acetate	M0631	Sigma-Aldrich	11.8	29.5	mM
Potassium glutamate	49601	Sigma-Aldrich	100	250	mM
DTT	sc-29089B	SantaCruz Biotech	1	2.5	mM
ATP	R0481	ThermoFisher	2	5	mM
GTP	R0481	ThermoFisher	2	5	mM
CTP	R0481	ThermoFisher	1	2.5	mM
UTP	R0481	ThermoFisher	1	2.5	mM
tRNA	10109541001	Roche	52	130	U _{A260} /mL
Creatine phosphate	27920	Sigma-Aldrich	20	50	mM
Folinic acid	F7876	Sigma-Aldrich	0.02	0.1	mM
Spermidine	S2626	Sigma-Aldrich	2	5	mM
HEPES	H0887-100ML	Sigma-Aldrich	50	125	mM